



50/10 GENESI M

Cod. D-GE150-040-M

Oil free compressor, FAD at 5 bar 120 l/min, tank capacity 40 L, with membrane air dryer and 0,01 micron filtration.


 Maximum simultaneous use: 2

 Dental Chairs: 2

Medical Device CE 0425

Air Dryer	Air capacity		FAD Fad @ 5bar		Motor	Voltage Hz	Receiver	Noise Level	Duty Cycle	max P.	Size	net Weight
Yes	m ³ /min	l/min	m ³ /min	l/min	kW	230V 50Hz	l	db(A)	50 %	bar	mm	Kg
	0.0064	6.4	0.12	120	1.1		40	65		7	710x410x770	42

Air delivery performance

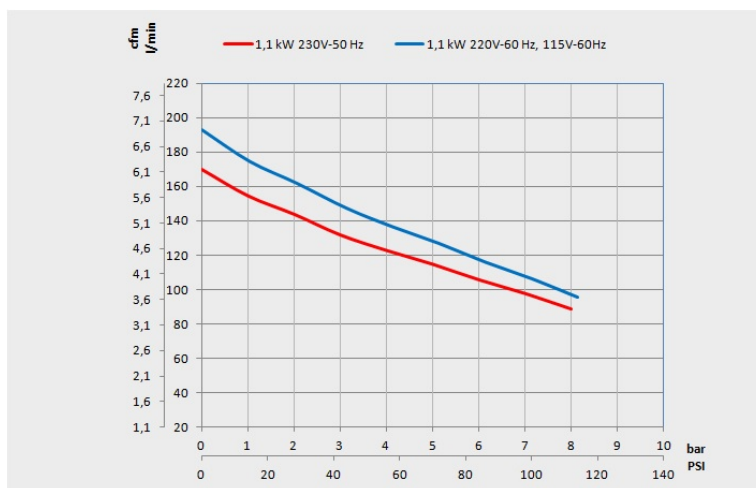

Maintenance kits

Kit.	Description	Type	Hours
TA-FL-024-K-SCC	<p>GENESI M ordinary maintenance kit</p> <p>Ordinary maintenance kit for GENESI compressors with dryer. The kit consists of the cartridges of the dryer filtration battery (2 cartridges: 5 + 0.01 micron) and the relative O-rings, the pumping element suction filter and the non-return valve pad. Annual replacement of kit components is recommended. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.</p>	Ordinary	1.000h/1y

Kit.	Description	Type	Hours
MG-GC-011-K	<p>Extra ordinary maintenance kit for GENESI S and M versions</p> <p>The GENESI extraordinary maintenance spare parts kit includes: a suction filter, two elastic rings, O-rings, a line filter cartridge and two complete valve plates. The recommended maintenance intervals for this compressor are 1500 hours / 2 years. Maintenance operations can be performed at a higher frequency if the environment in which the compressor works is particularly dirty.</p>	Extraordinary	1.500h/2y