

1st Brushless ENDO MOTOR

by Woodpecker

Ai-Motor

Super mini head & Slim neck
Stronger power
Precise rotation
Lower noise
Balances Weight




GUILIN WOODPECKER MEDICAL INSTRUMENT CO., LTD.




Scan and Login website for more information

High performance brushless motor

- 1- Adjustable reciprocating angle (10° interval; adjustment range: 20° -340°); Accurate rotating angle reduces risk of instrument fracture and enables higher preparing efficiency.

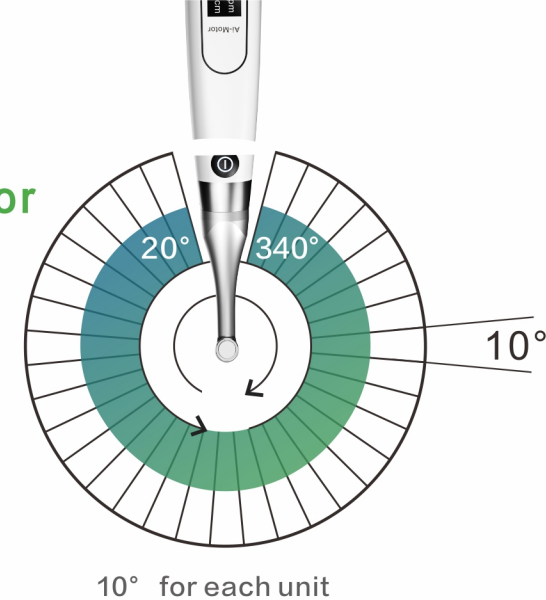

- 2- Smooth rotation and lower noise make treatment more endurable for patients.

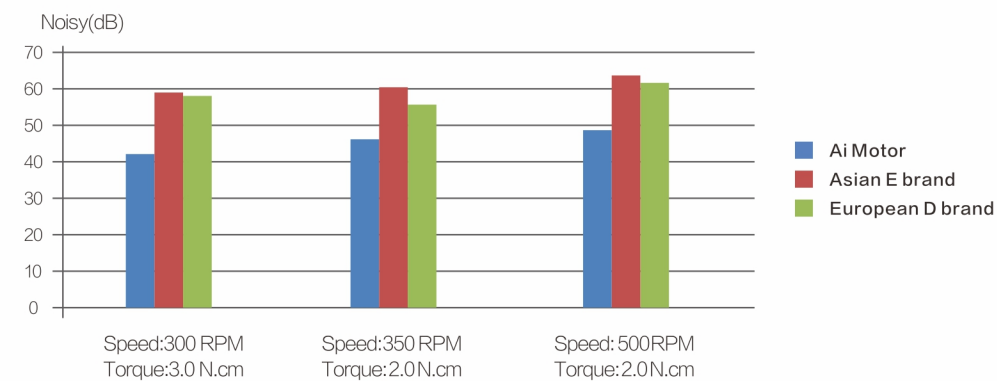


- 3- Strong power, together with accurate and stable torque output, helps to realize efficient root canal preparation.



360° rotatable contra-angle




2000mAh Lithium battery

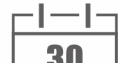

Built-in 30 mainstream file systems.

Super Mini Contra-angle


Super mini head with a diameter of 8mm and a height of 9.7mm brings clear vision during operation.

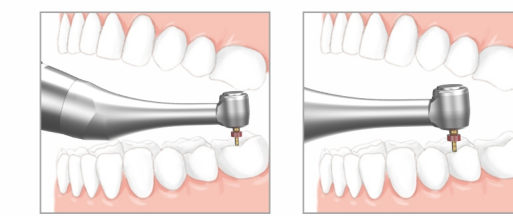

Woodpecker Super Mini head&Slim neck Normal head&neck

Built-in Apex locator

Integrated length determination makes endodontic treatment safer and more efficient. Realize real-time display of file location in canal. Automatically stops or reversely rotates as soon as the file reaches the apical stop, so as to prevent perforation. Auto Start and Auto Stop functions respectively enable motor to automatically start while the file gets in canal and stop when it is out.

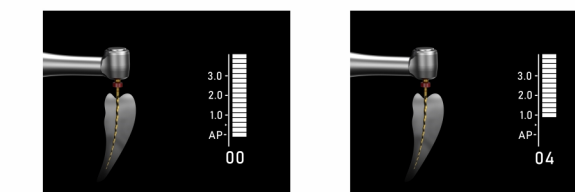

Convenient wireless charging



6:1 Contra-angle



Three times longer service life



Front center-of-gravity



The choice of millions of dentists worldwide

Guilin Woodpecker Medical Instrument Co., Ltd.
Information Industrial Park, Guilin National High-Tech Zone, Guilin, Guangxi, 541004 P. R. China

Tel:
Europe Sales Dept.: +86-773-5873196, +86-773-2125222
North America, South America & Oceania Sales Dept.: +86-773-5873198, +86-773-2125123
Asia & Africa Sales Dept.: +86-773-5855350, +86-773-2125896
Fax: +86-773-5822450
E-mail: woodpecker@glwoodpecker.com, sales@glwoodpecker.com
Website: http://www.glwoodpecker.com

EC REP MedNet GmbH
Borkstrasse 10 · 48163 Muenster · Germany