


**MY  
ITENA  
CHOICE**  
#SingleUse

**AIREO**

V-LOCK SYSTEM

**Single use tips for air/water syringes**

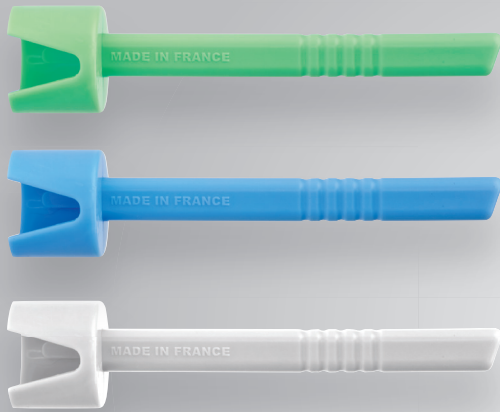

- **Guaranteed dry air:** air and water pass through two separate pipes right to the end of the tip
- **Single use:** avoids the risk of microbial migration
- **Reduces the risk of direct or cross contamination**
- «V» notch for easy insertion
- **Extremely flexible tip:** keeps the angle desired by the dentist
- Grooves improve flexibility
- Allows to reach hard-to-reach areas **without affecting the quality of the sprays**
- Made in France
- Recyclable medical grade plastic


### SETTING

- Simply push in the tip that fits perfectly to the adapter
- The nozzle holds firmly to the adapter to prevent the creation of a possible air chamber and therefore the loss of pressure at the outlet of the nozzle
- To change it, just pull

### References

200 AIREO V-Lock System tips box	
Green.....	AIRU200-V
Blue.....	AIRU200-B
White.....	AIRU200-A

### CIRCUIT LAYOUT

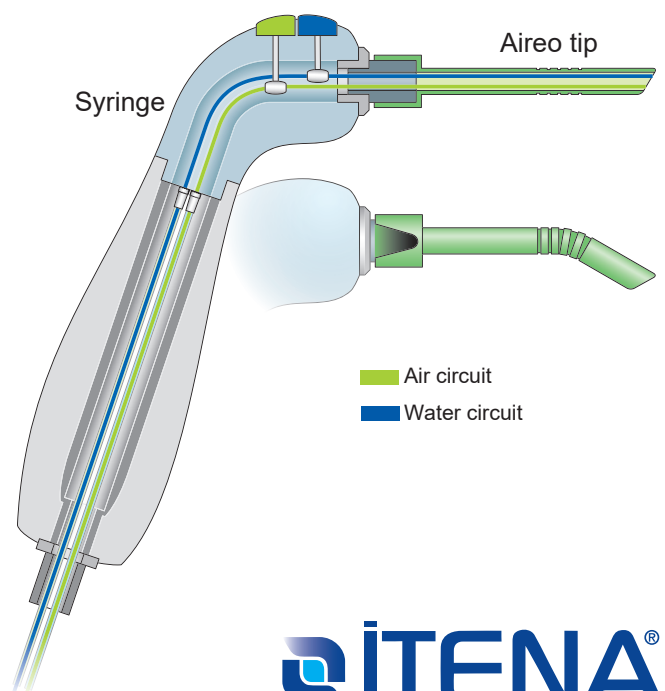

**ITENA**<sup>®</sup>  
CLINICAL PRODUCTS

Find us on [www.itena-clinical.com](http://www.itena-clinical.com)

Medical device marked CE. Reserved for dental professionals. Read instructions for use carefully before use.  
ITENA CLINICAL 31 avenue Georges Clemenceau, 93420 Villepinte - FRANCE - Tel: + 33 1 45 91 61 40 - IND OR - 11/2020